



# Evaluation of Coriolis Micro Air Sampling to Detect Volatile and Semi-Volatile Organic Compounds

Audrey Courtier, Benoit Roig, Stéphane Cariou, Axelle Cadiere, Sandrine Bayle

## ► To cite this version:

Audrey Courtier, Benoit Roig, Stéphane Cariou, Axelle Cadiere, Sandrine Bayle. Evaluation of Coriolis Micro Air Sampling to Detect Volatile and Semi-Volatile Organic Compounds. *Molecules*, 2022, 27 (19), pp.6462. 10.3390/molecules27196462 . hal-03821482

**HAL Id: hal-03821482**

**<https://imt-mines-ales.hal.science/hal-03821482>**

Submitted on 19 Oct 2022

**HAL** is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

## Article

# Evaluation of Coriolis Micro Air Sampling to Detect Volatile and Semi-Volatile Organic Compounds

Audrey Courtier <sup>1,\*</sup> , Benoit Roig <sup>1</sup>, Stephane Cariou <sup>2</sup>, Axelle Cadere <sup>1</sup>  and Sandrine Bayle <sup>2</sup> 

<sup>1</sup> UPR Chrome, University of Nimes, Rue du Dr G. Salan, CEDEX 1, 30021 Nimes, France

<sup>2</sup> Laboratoire des Sciences des Risques (LSR), IMT Alès, 6 Av. de Clavières, 30100 Alès, France

\* Correspondence: audrey.courtier@unimes.fr

**Abstract:** There are several analytical procedures available for the monitoring of volatile organic compounds (VOCs) in the air, which differ mainly on sampling procedures. The Coriolis micro air sampler is a tool normally designed for biological air sampling. In this paper, the Coriolis micro bio collector is used to evaluate its ability to sample organic contaminants sampling and detecting them when combined GC-MS. We also compare the use of the Coriolis micro with a standardized sampling method, which is the use of a lung box with a Nalophan<sup>®</sup> bag. The results show that the Coriolis micro sampling method is suitable for the sampling of organic contaminants. Indeed, the Coriolis micro allows to sample and detect mainly semi-volatile molecules, while the lung box/Nalophan<sup>®</sup> bags allow to sample more volatile molecules (highly volatile and volatile). These results were confirmed in the controlled air lab with a slight difference with the field. The simultaneous use of the both techniques allow to sample and detect a larger number of molecules with specific physicochemical properties to each sampling technique. In conclusion, the Coriolis micro can sample and detect volatile organic compounds present in air. We have shown that the development of alternative sampling methods and the use of non-target analysis are essential for a more comprehensive risk assessment. Moreover, the use of the Coriolis micro allows the detection of emergent molecules around the Thau lagoon.

**Keywords:** atmospheric pollutants; Coriolis micro; Nalophan<sup>®</sup> bags; environmental detection



**Citation:** Courtier, A.; Roig, B.; Cariou, S.; Cadere, A.; Bayle, S. Evaluation of Coriolis Micro Air Sampling to Detect Volatile and Semi-Volatile Organic Compounds. *Molecules* **2022**, *27*, 6462. <https://doi.org/10.3390/molecules27196462>

Academic Editor:  
Carmen González-Barreiro

Received: 3 August 2022  
Accepted: 22 September 2022  
Published: 30 September 2022

**Publisher's Note:** MDPI stays neutral with regard to jurisdictional claims in published maps and institutional affiliations.



**Copyright:** © 2022 by the authors. Licensee MDPI, Basel, Switzerland. This article is an open access article distributed under the terms and conditions of the Creative Commons Attribution (CC BY) license (<https://creativecommons.org/licenses/by/4.0/>).

## 1. Introduction

The increase of industrialization and urbanization leads to growing pollution in different compartments such as air, water, and soil. These contaminations can induce risks for human health and the environment. The World Health Organization (WHO) alerts that 13.7 million deaths per year are related to the environment and that 4.2 million of these are the result of exposure to outdoor air pollution [1]. Outdoor air pollution has been associated with many diseases such as childhood asthma [2], kidney diseases [3], lung cancer [4], and cardiopulmonary disease [5].

The most common outdoor and indoor air pollutants are volatile organic compounds (VOCs) which are defined as any compound of carbon in atmospheric photochemical reactions, according to the United States Environmental Protection Agency (US-EPA) [6]. VOCs can be natural and anthropogenic sources and some of them are known to be of health concern [7–11]. Outdoor air pollutants are a heterogeneous class of molecules that include particles such as bioaerosols, semi volatile compounds (SVOCs), volatile (VOCs), and very volatile compounds (VVOCs). Categories of chemical pollutants are defined by specific properties. For example, a VOC is defined by a boiling point (BP) ranging from 50–100 °C to 240–260 °C, while a SVOC has a BP ranging from 240–260 °C to 380–400 °C according to WHO [12]. The WHO intentionally does not clearly define between VVOCs and VOCs, while the US-EPA defined VVOCs as compounds having a BP lower than 50–100 °C [13]. Particles are quantified via airborne particulate matter (PM) and may represent a danger to

human health. PM consists of a heterogeneous mixture of airborne solid and liquid particles that varies continuously in size and chemical composition. The chemical constituents of PM are highly heterogeneous and include nitrate, sulfates, elemental and organic carbon, organic compounds, biological compounds, and metals [14].

The greatest challenge in humans and/or the environment assessing risk is to sample as many airborne pollutants as possible that may pose a risk. Several analytical methods for the determination of VOCs in the air exist, which differ mainly by the sampling procedures [15]. The sampling procedures consist of passive or diffusive sampling, active sampling, the use of canisters and bags, and online sampling [16]. Passive sampling is based on the molecular diffusion of molecules through a surface and does not detect the pollution peak. It allows the sampling of most VOCs but can also sample SVOCs and particles [17,18]. In some studies, passive sampling has been used for VVOCs but has shown an underestimation of their concentration [19]. Active sampling is preferable for sampling VVOCs [20]. Active sampling consists of passing a constant flow of air through a tube filled with sorbent to catch VOCs. Sorbents can be of different types (activated carbon, polymers, silica, etc.) with one or several types of sorbents [21]. The main difficulty of this type of sampling is to maintain the highly volatile molecules trapped on the sorbents and the inherent artifact, depending on the sorbents used [21]. A canister consists of an inert stainless steel spherical container of several liters (1 to 15 L). After sampling, an aliquot of gas is transferred to a pre-concentrator trap of thermal desorber. It is also possible to collect air in bags and different types are available, such as Tedlar<sup>®</sup>, Kynar<sup>®</sup> or Nalophan<sup>®</sup> bags. Canisters or bags are a useful alternative to sorbent tubes, particularly for ultra-volatile compounds [16], and they can also be used for sampling particles [22]. The last method, that is the one currently used, is online sampling. The air samples are drawn manifold directly into the sorbent focusing trap of the thermal desorber. This technique requires a complete analytical system at each sampling point, but reduces the errors due to possible reactions occurring during storage of the adsorbent tubes before analysis [16].


Each sampling method has its advantages and disadvantages, as discussed above. The challenge lies especially in the sampling of SVOCs and VVOCs, i.e., managing to sample a large number of molecules with different properties. After sampling, two types of analysis can be performed to identify contamination: targeted and/or non-targeted analyses. Targeted analysis is generally used in order to find specific pollutants, which are mostly regulated substances [23,24]. Conversely, non-targeted analyses allow the discovery of new contaminants of emerging concern without “a priori”. Currently, environmental health studies, particularly for risk assessment, emphasize non-targeted approaches. These techniques allow the detection of many molecules [25–27]. Consequently, the results provided by these techniques highlighted new molecules whose health impact was previously unsuspected [28]. The right combination of sampling and analytical method allows a more accurate risk analysis.

In order to detect a maximum of compounds around the study site, a new methodology, including an original sampling/analysis combination, was developed. For the sampling method, the Coriolis micro (Bertin instruments, France), normally used to concentrate microorganisms in the air, was assessed. Indeed, this device has the advantage of concentrating air particles in water at a high flow rate, reducing the sampling time. Moreover, this technique allows the recovery of hydrophilic molecules or SVOCs. To extract and concentrate compounds from the water, a SBSE/HSSE extraction was realized before performing a non-targeted analysis by GC/MS (Agilent 7890B gas chromatograph/Agilent 5977A mass spectrometer). These analyses are compared to a standardized sampling method (NF X 43-104 from 1995 recently adopted by the EN 13725:2022 standard) which is a lung box containing Nalophan<sup>®</sup> bags, followed by desorption on multi-bed sorbent tubes, and a non-targeted GC-MS analysis. Comparisons were performed in an operational context, either in the field or under controlled conditions.

## 2. Results

### 2.1. Comparison of the Two Sampling Techniques

The air around the Thau lagoon was collected for 40 min using the lung box/Nalophan<sup>®</sup> bag and the Coriolis micro techniques. The volumes of air collected were 40 L and 12,000 L, respectively. The air volume was higher with the use of the Coriolis. Fifty-two molecules were found with the Nalophan<sup>®</sup> bags, while forty-four were found with Coriolis micro. Two molecules, nonanal and 2-ethyl-1-hexanol, were found in common (Figure 1). The volume of air sampled does not impact the number of molecules obtained but rather the nature of the molecules. Furthermore, it seems that the combination of both techniques allows to obtain a greater number of molecules. To understand which types of molecules were sampled according to the technique used, the physico-chemical properties of the molecules recovered with both techniques were sought (Supplementary Materials: Tables S1 and S2).

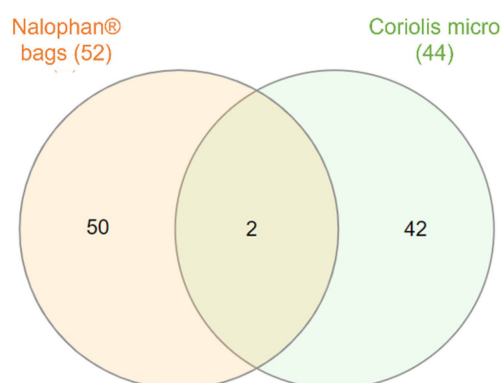

**Figure 1.** Distribution of molecules found using both techniques of air sampling. Venn diagram illustrating the degree of overlap of molecules found with Coriolis micro and Nalophan<sup>®</sup> bags method. The central section in orange represents the compounds found in common between the two methods. Specific compounds found with Nalophan<sup>®</sup> bags are in pink, whereas those found with Coriolis micro are in blue.

### 2.2. Compounds Properties

To assess the types of molecules retained by each of the two techniques, different physico-chemical parameters were investigated: molecular weight, the behavior in water compartment with solubility at 25 °C and the octanol/water partition coefficient, the behavior in air compartment with the octanol/air partition coefficient, Henry's constant, boiling point, and vapor pressure.

#### 2.2.1. Molecular Weight


The comparison of the molecular weight of the molecules detected by the two techniques is significantly different for each sampling point (all  $p$ -Value < 0.01, see Table S3). The weight is higher for the Coriolis sampling method with a mean at 156 g/mol than the Nalophan<sup>®</sup> bag, with a mean at 103 g/mol (Figure 2).

#### 2.2.2. Behavior in Water Compartment

- Solubility in Water at 25 °C

The solubilities of the molecules found with each technique are different. Indeed, the solubility range is greater with the Coriolis micro than with the Nalophan<sup>®</sup> bag, with no significant difference except for P6, while the medians are quite similar (Figure 3a). It is interesting to note the solubilities of the molecules obtained with the Coriolis micro are not higher than those obtained with the Nalophan<sup>®</sup> bags, which could have been the case since the molecules sampled with the Coriolis micro are recovered in water.

- The Octanol/Water Partition Coefficient ( $\log K_{ow}$ )

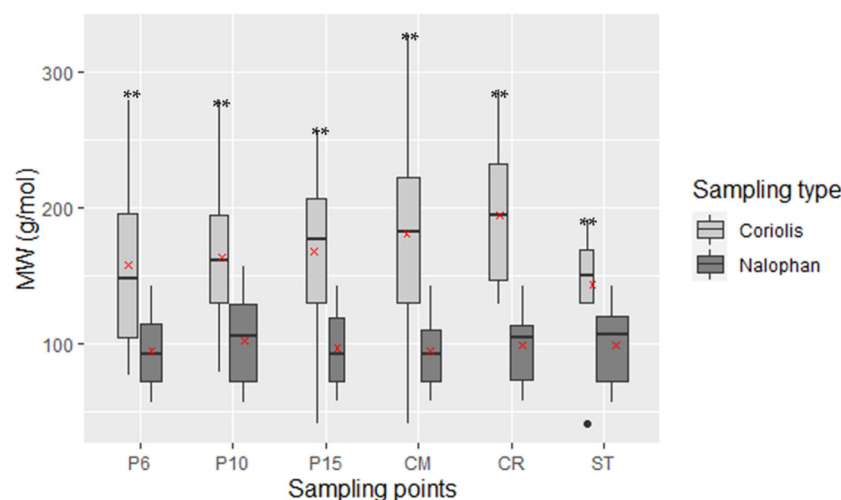

**Figure 2.** Comparison of molecular weights according to the sampler type used. The red crosses represent the mean values. MW are significantly different between two techniques (\*\* mean  $p$  value < 0.01) values follow normal distribution except ST values. Black dott represents an outlier value of the boxplot.

To better understand the capacity to these molecules to bioaccumulate in organisms, the water/octanol partition coefficient ( $\log K_{ow}$ ) was inquired.  $\log K_{ow}$  determines if molecules are more hydrophobic or lipophilic. In other words,  $\log K_{ow}$  reflects the distribution capacity of organic compounds between octanol (hydrophobic/lipophilic) and water (hydrophilic). Molecules with  $\log K_{ow}$  less than 3 are considered hydrophilic compounds, whereas molecules with  $\log K_{ow}$  higher than 3 are considered hydrophobic compounds.

The boxplot shows a wide range of  $\log K_{ow}$  from  $-2.3$  to  $7$  for Coriolis micro and from  $-1.2$  to  $4.8$  for Nalophan<sup>®</sup> bags, the widest being for the Coriolis micro but with no significant difference (Figure 3b). However, the medians are almost identical with  $2.69$  for the Coriolis micro and  $2.65$  for Nalophan<sup>®</sup> bags. The same proportion of soluble and/or lipophile molecules was obtained with  $66\%$  for the Coriolis micro and  $60\%$  with Nalophan<sup>®</sup> bags.

### 2.2.3. Behavior in Air Compartment


- The Octanol/Air Partition Coefficient ( $\log K_{oa}$ )

The  $\log K_{oa}$  is useful to predict the partitioning behavior of organic compounds between air and environmental matrices, such as soil, vegetation, and aerosol particles.  $\log K_{oa}$  is the ratio of a chemical's concentration in octanol to the concentration in air equilibrium. Organic molecules with  $\log K_{oa}$  comprised between  $6.5$  and  $10$  have a significant potential for long-distance air transport, whereas organic compounds with  $\log K_{oa} < 6.5$  have a high volatility and cannot be deposited on the ecosystem. Finally, organic compounds with  $\log K_{oa} > 10$  have a greater potential for adsorption by the particles surface and lipids [29,30].  $\log K_{oa}$  also gives information about bioaccumulation. Molecules with  $6 < \log K_{oa} < 12$  have a high bioaccumulation potential for humans [30].

Figure 4a highlights  $\log K_{oa}$  repartition for each method and shows repartitions ranging from  $2.494$  to  $10.603$  for the Coriolis micro and from  $0.712$  to  $6.255$  for Nalophan<sup>®</sup> bags. The medians are different with  $6.26$  and  $3.44$  for the Coriolis micro and Nalophan<sup>®</sup> bags, respectively. Values obtained for Nalophan<sup>®</sup> bags are significantly different from the Coriolis values ( $p$ -values < 0.01).  $\log K_{oa} < 6.5$ , which defined high volatility, characterized all the molecules obtained with Nalophan<sup>®</sup> bags against only  $52\%$  of molecules obtained with the Coriolis micro. Molecules sampled with the Coriolis micro had other particularities, namely potential for long distance air transport ( $6.5 < \log K_{oa} < 10$ ), which concern  $43\%$  of the molecules. A small proportion of molecules obtained with the Coriolis micro ( $4.5\%$ )

had  $\log K_{oa} > 10$ , which means a greater potential to be adsorbed by the particles surface and lipids.

In conclusion, molecules found with the Coriolis micro had higher  $\log K_{oa}$  than molecules identified with Nalophan® bags. This means that molecules found with the Coriolis had potential for long-distance air transport, bioaccumulation in humans, and adsorption to air particles. Molecules found with Nalophan® bags were more volatile than molecules found with the Coriolis.

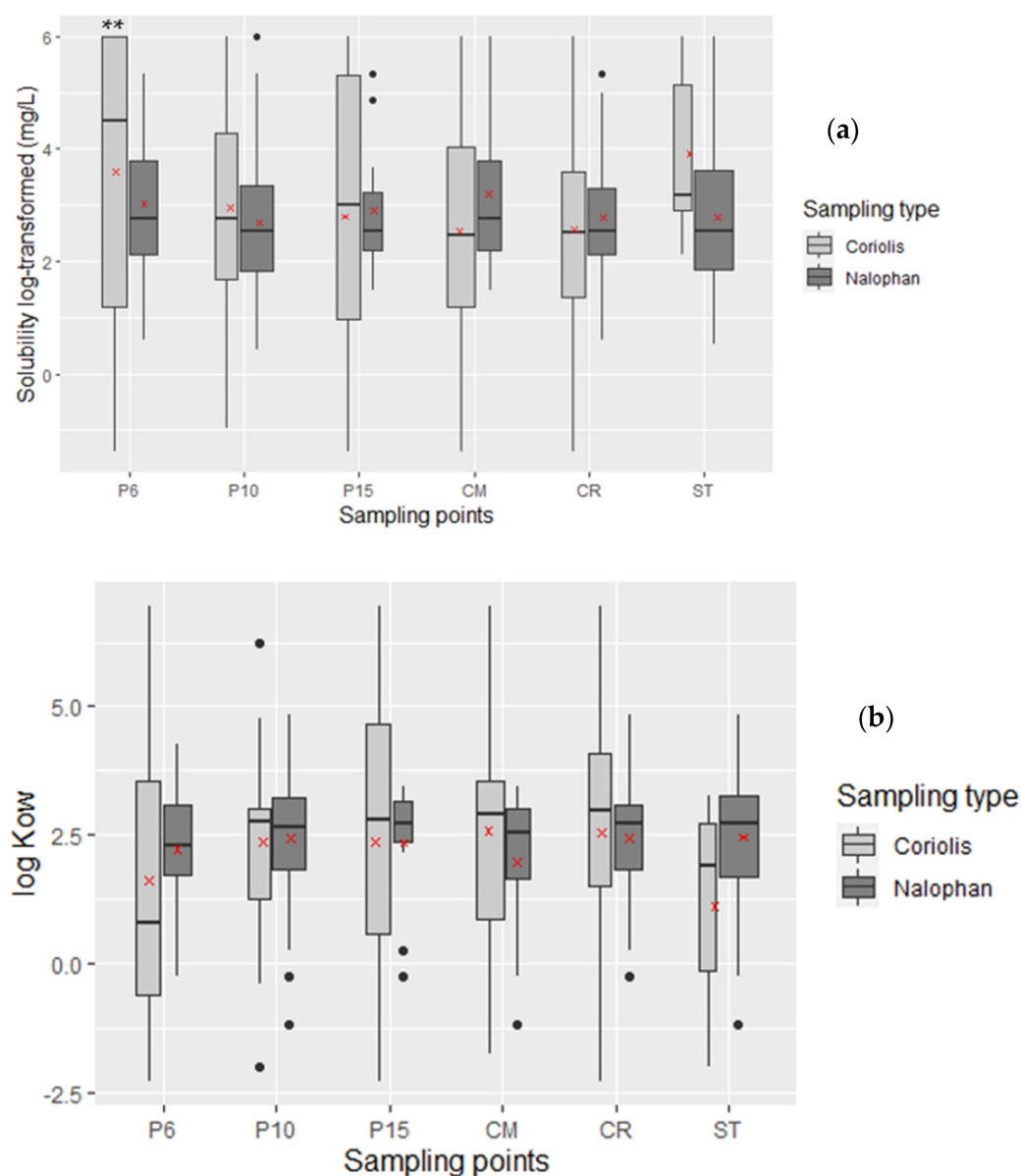

**Figure 3.** Comparison of: (a) solubilities in water at 25 °C according to the sampler type used. The red crosses represent the mean values. No significant difference is found between these two variables except for P6; (b) the octanol/water partition coefficient according to the method used. No significant difference is found between these two variables. \*\* mean  $p$  value < 0.01. Black dott represents an outlier value of the boxplot.

- Henry's Constant

Henry's constant characterizes the capacity of a substance to split between the two phases of a binary air/water system, Henry's constants higher than  $10 \text{ Pa}\cdot\text{m}^3/\text{mol}$  allows to predict a volatilization of the substance from water bodies, while Henry's constant com-

prised between 0.1 to 10 Pa.m<sup>3</sup>/mol allows to predict a slow volatilization [31]. Molecules retained by the two techniques are represented in function of their Henry's constant distribution. Box plot highlighted that the individual's repartition is not significantly different for each method, except for P6 ( $p$ -value < 0.01), even if Nalophan<sup>®</sup> bags values are higher than values obtained for the Coriolis (Figure 4b). It reveals that median values were different with 2.9 and −0.53 for Nalophan<sup>®</sup> bags and the Coriolis micro, respectively. The Henry's constant repartition shows that 84% of molecules found with the Coriolis micro appears to have a Henry's constant lower than 10 Pa.m<sup>3</sup>/mol (equal to 1 in log-transformed on the figure), which means slow volatilization, against 21% with Nalophan<sup>®</sup> bags (data not shown). Conversely, 16% and 79% of molecules found with the Coriolis micro and Nalophan<sup>®</sup> bags, respectively, had a Henry's constant greater than 10 Pa.m<sup>3</sup>/mol. Consequently, molecules found with Nalophan<sup>®</sup> bags are more volatile than molecules found with the Coriolis micro.


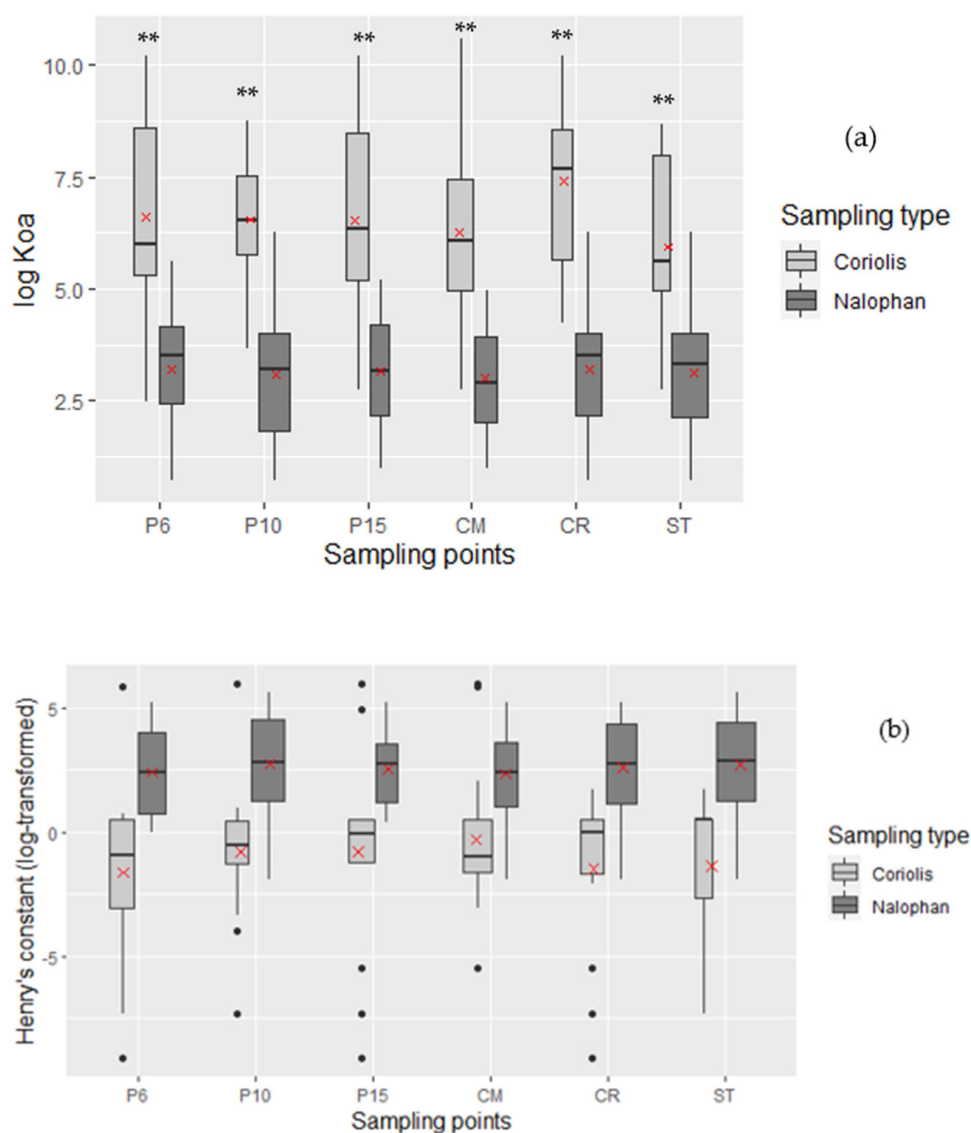


Figure 4. Cont.

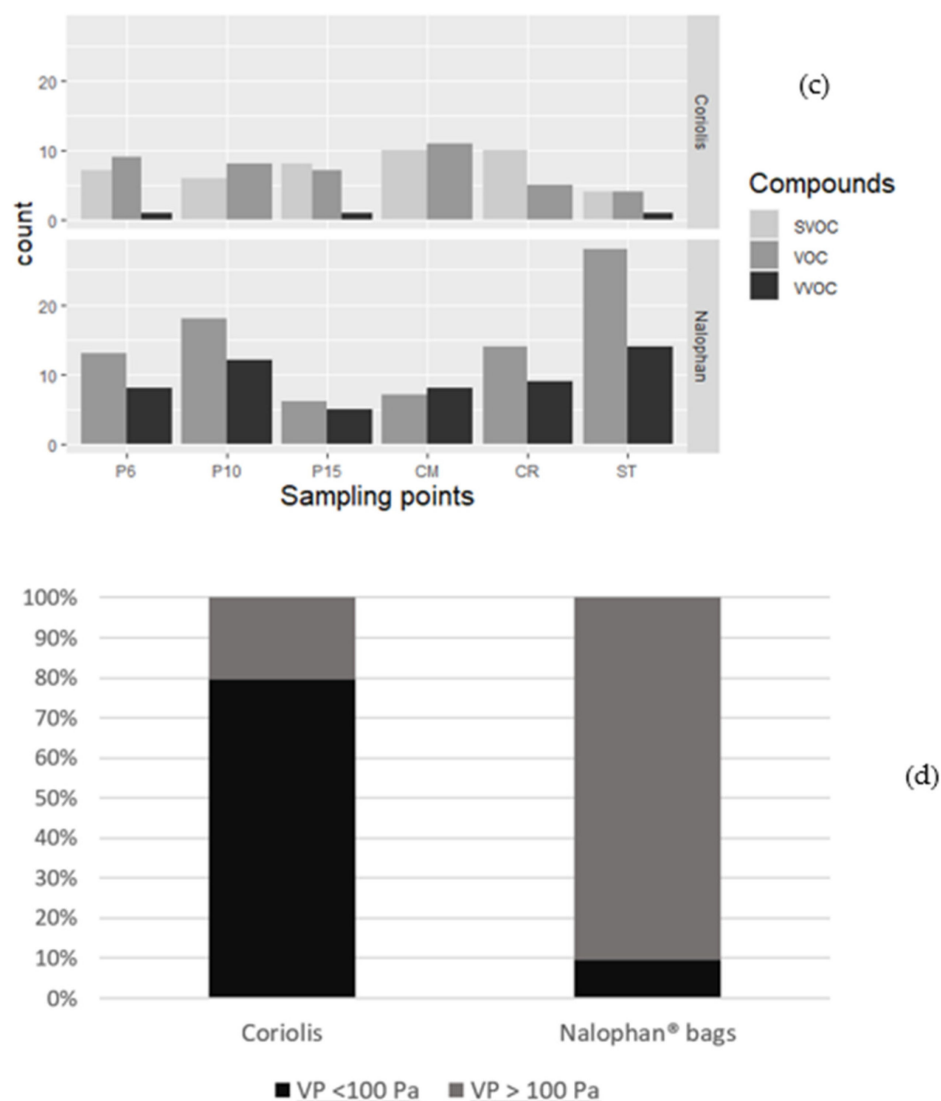

**Figure 4.** (a) The octanol/air partition coefficient repartition according to the method used. Log  $K_{oa}$  of molecules obtained with Coriolis micro are significantly different from those obtained with Nalophan® bags ( $p$ -value < 0.01, Supplementary Materials: Table S3); (b) box plot of the log-transformed Henry's constant in  $\text{Pa}\cdot\text{m}^3/\text{mol}$  molecules according to the methods used. Henry's constants are not significantly different between those of Coriolis micro and those of Nalophan® bags except for P6 ( $p$ -value < 0.01); (c) repartition of boiling points (BP) according to the method use. Figure represents BP comprised between 260 °C and 400 °C for semi volatile organic compound (SVOC); BP comprised between 100 °C and 260 °C for volatile organic compound (VOC) and BP below 100 °C for very volatile organic compound (VVOC), no significant differences were found; (d) vapor pressure superior and inferior to 100 Pa to determine volatility of each compound found with both techniques, no significant differences were found. Black dott represents an outlier value of the boxplot; \*\* mean  $p$  value < 0.01.

- Boiling Point (BP)

A VOC is defined with a BP ranging from 50–100 °C to 240–260 °C, according to WHO [12], and a semi volatile organic compound (SVOC) has a BP ranging from 240–260 °C to 380–400 °C. A very volatile organic compound (VVOC) has a BP lower than 50–100 °C [13].

All molecules found with Nalophan® bags had a BP below 260 °C (Figure 4c). Indeed, 35 molecules found with Nalophan® bags had BP comprised between 100 °C and 260 °C,

which is a characteristic of VOC, while 17 molecules had a BP below 100 °C, which corresponds to VVOC. On the other hand, 21 molecules found with the Coriolis micro had a BP comprised between 260 °C and 400 °C, which is a characteristic of semi volatile molecules. However, no significant differences were obtained between values of Nalophan<sup>®</sup> bags and the Coriolis micro. Moreover, the same number of molecules had a BP comprised between 100 °C and 260 °C and only 2 with a BP lower than 100 °C. In other words, the Coriolis micro enables to sample SVOC, VOC, and slightly VVOC, while Nalophan<sup>®</sup> bags permit to samples mainly VOC and VVOC.

- Vapor Pressure (VP)

A vapor pressure greater than 100 Pa characterizes a volatile compound and a vapor pressure comprised between 0.001 and 100 Pa characterizes a semi volatile compound [32].

There are 79.6% of the Coriolis micro molecules with vapor pressure <100 Pa, against 10% for Nalophan<sup>®</sup> bags. Inversely, 90.4% of Nalophan<sup>®</sup> bags molecules had a vapor pressure >100 Pa (Figure 4d). This supports the fact that the molecules obtained with the Nalophan<sup>®</sup> bags are considered as volatile, while the Coriolis micro allows to sample semi-volatile molecules.

These parameters allow to confirm that molecules found with Nalophan<sup>®</sup> bags are more volatile than the ones found with the Coriolis. Molecules found with the Coriolis are defined as semi volatile compounds and can have a more complex behavior. On the other hand, molecules obtained with the Coriolis micro can be adsorbed on particles and potentially fall back into a compartment different from their origin. These molecules can remain close to their initial site or be subjected to long distance air transport.

### 2.3. Sampling in Controlled Air


To better understand the type of molecules collected by each of the two methods, we carried out samples in controlled air. These samplings allow us to avoid the variations which can be present on the field, and can be different according to the sampling points (wind, temperature, etc.)

A sampling test with the Coriolis micro and Nalophan<sup>®</sup> bags was carried out in controlled indoor air. Seven molecules were present in a standardized solution used to make artificial contaminated air. Among these molecules, two were found in the field only with the Coriolis micro (1-octanol and undecane), two were found in the field only by Nalophan<sup>®</sup> bags (toluene and 2-butanone), and two were found by both techniques (octanal and 2-ethyl-1-hexanol). The results show that the Coriolis micro sampled three molecules out of the four sampled in the field (Table 1). Among the expected molecules, only undecane was not sampled with the Coriolis micro in controlled air. Indeed, the Coriolis micro has sampled 1-octanol, octanal, and 2-ethyl-1-hexanol, in all samples ( $n = 3$ ). Nalophan<sup>®</sup> bags has sampled all the molecules present in the solution (7 out of 7 molecules). 1-octanol and undecane were not sampled by Nalophan<sup>®</sup> bags in the field, while they were detected with the Coriolis micro. Sampling in controlled air showed a slight lack of sensitivity for the Coriolis micro, which did not detect a molecule present in the mixture in the three samples. Concerning Nalophan<sup>®</sup> bags, all the molecules were found in controlled air but not in the field. This may imply a higher detection limit not allowing to detect some compounds presents in the environment, or a matrix effect. The selected molecules have their physico-chemical properties in the range of the detection capacities of each technique, as seen previously (Figure 5). All the properties among partage coefficient octanol/air, octanol/water, Henry's constant, molecular weight, boiling point, vapor pressure, and solubility, seem to influence the ability of the Coriolis micro or Nalophan<sup>®</sup> bags to sample a type of molecule. Indeed, all the molecules are included in the box plot, except one extreme, which is undecane for both techniques and toluene and 2-butanone for the Coriolis micro. Undecane has a low solubility in water (0.03 mg/L), which can explain why we do not find it using the Coriolis micro. The undecane was already found in the air using the Coriolis micro, but only twice in the same campaign (9 other campaigns were developed within the framework of another study and no undecane was detected: data not shown).

Octanal is a widely present molecule that can be found in high concentrations in limes, mandarins, oranges, and other different foods [33,34]. Aldehyde (and octanal) are also known to be uptaken by acidic particles and form secondary organic aerosols (SOA), but this phenomenon implies a high concentration of octanal present in the atmosphere [35,36]. Moreover, it is important to note that the controlled air tests do not simulate the presence of particles in the air, while the Coriolis micro is used to sample the (bio)particles.

**Table 1.** Molecules detected (yes) or not (no) according to sampling type and place (lab/field), ( $n = 3$ ).

| Molecules         | Field          |                | Laboratory     |                |
|-------------------|----------------|----------------|----------------|----------------|
|                   | Coriolis Micro | Nalophan® Bags | Coriolis Micro | Nalophan® Bags |
| 1-octanol         | yes            | no             | yes            | yes            |
| 2-butanone        | no             | yes            | no             | yes            |
| 2-ethyl-1-hexanol | yes            | yes            | yes            | yes            |
| Octanal           | no             | yes            | yes            | yes            |
| Toluene           | no             | yes            | no             | yes            |
| Undecane          | yes            | no             | no             | yes            |

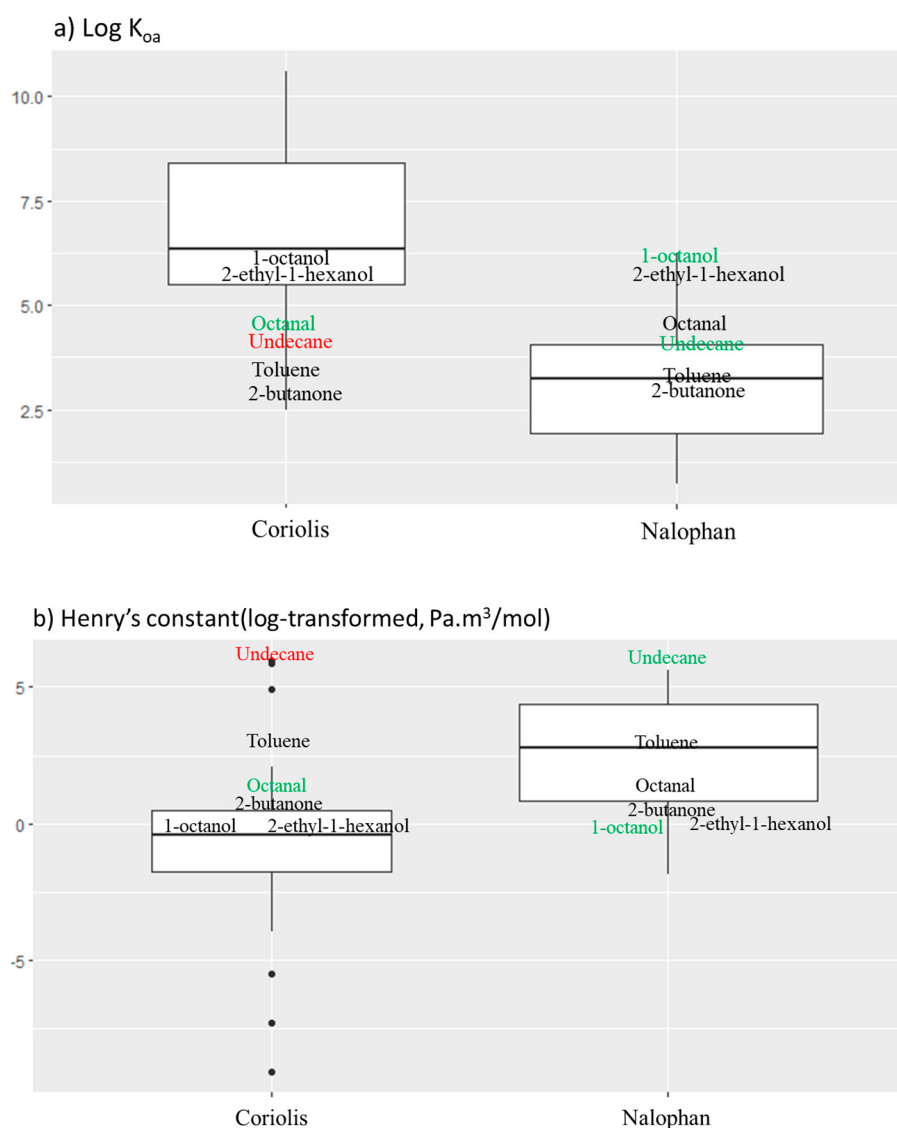

**Figure 5.** Cont.


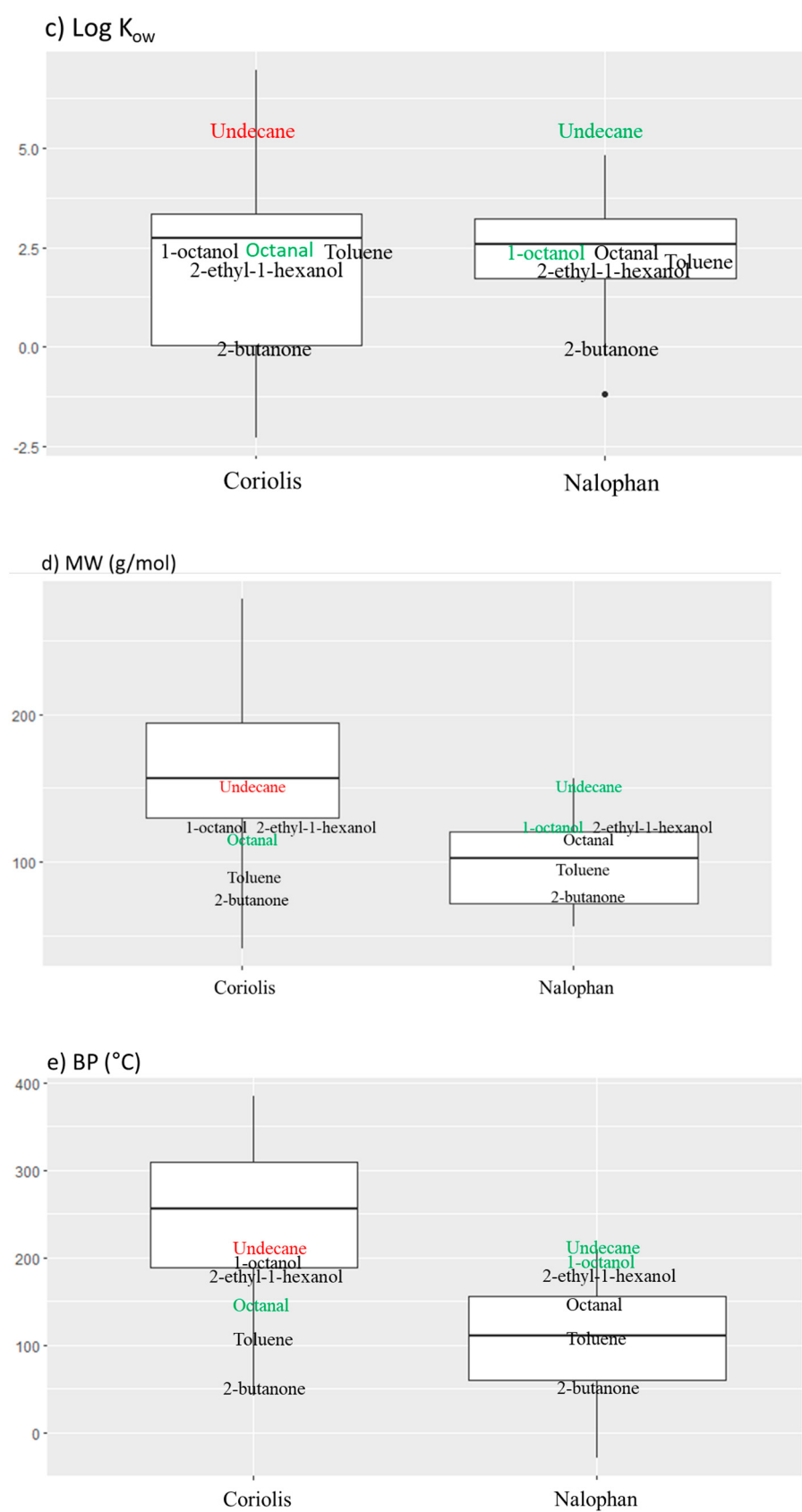


Figure 5. Cont.

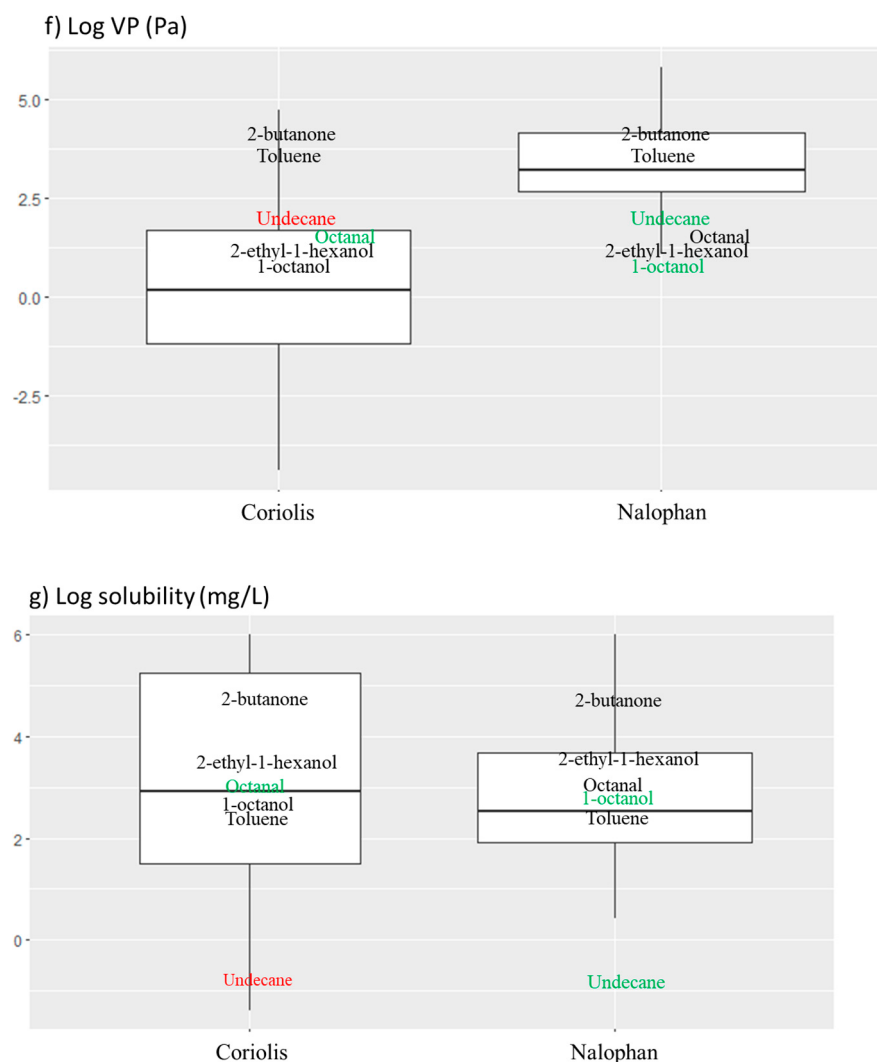

**Figure 5.** Box plot for each property previously presented: (a) log  $K_{oa}$ , (b) log-transformed Henry's constant, (c) log  $K_{ow}$ , (d) molecular weight (MW), (e) boiling point (BP), (f) log vapor pressure (VP), and (g) log-transformed solubility (mg/L) with position of each molecule present in solution used to make artificial contaminated air. The molecules in green color are the molecules found only during the laboratory test, the molecules in red color are the molecules found only during the field. Black dotted represents an outlier value of the boxplot.

### 3. Discussion

A sampling method usually dedicated to the collection of microorganisms was tested. The technique is a cyclonic impaction method. According to the manufacturer's information, particles are collected from 0.5  $\mu\text{m}$ . The compounds detected by this method can be either compounds dissolved into the water or compounds bound to the particles. It was expected that the evaluated method would retain the polar and less volatile molecules, with possibly other compounds retained on the particles, while the reference method would present the more volatile molecules. However, it is not easy to provide the impact of sampling a large volume of air. The manipulations proposed in our study allow us to confirm our hypothesis.

The results of this study showed a great complementarity in the two sample methods in terms of physico-chemical properties, such as molecular weight and the octanol/air partition coefficient. More specifically, the Coriolis micro allows to sample SVOCs, VOCs, and to a lesser extent VVOCs. The molecules obtained with the Coriolis micro have shown other properties, such as long-distance potential and adsorption on particles. Additionally, Nalophan<sup>®</sup> bags were able to collect VOCs and VVOCs.

Among the molecules found with both techniques, some of them are classified as of concern and require further investigation to assess their possible impact on the environment.

Among the emerging pollutants of interest, two in particular have attracted our attention, such as pyridine or methyl vinyl ketone (MVK). Pyridine is ubiquitous and is used as a solvent and intermediate in the production of agricultural chemicals, drugs, and paints [37]. Pyridine is also used in other processes i.e., a derivatization agent associated with methyl imidazole for converting non-volatile substances in volatile targets [38]. Exposure to pyridine and its derivatives, such as 3-picoline, 2-aminopyridine, 4-aminopyridine, etc., has harmful effects on the liver, kidneys, reproductive functions, and the immune system. Furthermore, it has potential carcinogenicity and has toxicity to aquatic life [39]. Pyridine has a very high half-life when released in the atmosphere [40] and it is also included in the list of chemicals subjected to EPCRA in the section, toxic chemicals [41]. MVK is one of the two major compounds produced during the oxidation of isoprene by the OH and is ozone-initiated, which is the largest atmospheric emission flux among all non-methane VOCs [42]; it may also have anthropogenic sources [43]. MVK can be hazardous to human health with upper respiratory tract effects [44] and is on the consolidated list of chemicals subject to the emergency planning and community right-to-know act (EPCRA), in the section, extremely hazardous substances, from USEPA [44].

Currently, atmospheric exposure risk assessments are performed on a limited number of methods, most often using targeted analyses [45–49]. Non targeted analyses are beginning to be used, particularly for the food matrices [50,51], and have shown their usefulness in facilitating provisional safety assessments and in finding new emerging chemicals [28,52]. Besides the non-targeted analysis, we show here that the development of alternative sampling methods is also essential for a more comprehensive risk assessment. The use of the Coriolis has allowed to highlight other molecules of interest and those emerging, such as phthalates (benzyl butyl phthalate, diisobutyl phthalate, dimethyl phthalate, and diethyl phthalate), phenol, pyridine, triethyl phosphate, and 1,2-benzothiazole, among other.

## 4. Materials and Methods

### 4.1. Sampling Points

The sampling was carried out in the south of France around the Thau lagoon. As shown in Figure 6, six different sampling points were selected around the lagoon, according to other studies previously developed in the framework of regulatory monitoring [53,54]. Among several sampling sites (from P1 to P20), three points which are located on the main rivers entering the lagoon were selected. P6 corresponded to the end of the Fontanilles stream, P10 corresponded to the end of the Pallas stream, and P15 corresponded to the end of the Vène stream. To continue in the follow-up of the major tributaries arriving in the lagoon, we decided to add the point CM, corresponding to the end of the canal of midi, and the point CR, corresponding to the beginning of the canal of Rhône. Finally, a last sampling point, directly in the lagoon and near the tourist town of Sète, was chosen and called ST.

### 4.2. Samples Collection and Analyses

At each point, air samples were collected using two different techniques. The first one consists of collecting air using laboratory-made 40 L Nalophan® bags and a lung chamber, according to the NF X 43-104 air sampling standard. The second technique involves using the Coriolis Micro (Bertin instruments, France) to collect air particles in 15 mL of distilled water for 40 min at 300 L/min. These sampling conditions correspond to the usual conditions of use in the field. The result is a difference in the sampling volume. Then, extractions and analyses were performed on each type of sample (Table 2). One sample was sampled at each site by each technique in a duplicate experiment.



**Figure 6.** Distribution of sampling points around Thau lagoon.

**Table 2.** Methods collection and analysis used for air samples.

| Stages     | Description                | Nalophan® Bag                                                                                                 | Coriolis Micro                                                                                                                                                                                                                                                                                                  |
|------------|----------------------------|---------------------------------------------------------------------------------------------------------------|-----------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Sampling   | Instruments                | Laboratory-made 40 L Nalophan® bag + lung chamber (NF X 43-104)                                               | Coriolis micro (Bertin instruments, France)                                                                                                                                                                                                                                                                     |
|            | Quantity of air sampled    | 40 L                                                                                                          | 12,000 L                                                                                                                                                                                                                                                                                                        |
| Extraction | Recovery unit              | 18 L are pumped through AirToxic (Perkin Elmer, Waltham, MA, USA) multi-bed sorbent tubes                     | 15 mL of distilled water in an amber glass vial. SBSE <sup>a</sup> with 10 mm Twister (Gerstel, Germany) magnetic stir bars coated with polydimethylsiloxane (PDMS) (0.5 mm film thickness) or ethylene glycol Silicone (EGS) <sup>*</sup> . The method was an adaptation of the Berrou et al., 2020 study [55] |
|            | Flow rate                  | 100 mL/min                                                                                                    | -                                                                                                                                                                                                                                                                                                               |
|            | Agitation-duration         | -                                                                                                             | 600 rpm during 2 h                                                                                                                                                                                                                                                                                              |
|            |                            |                                                                                                               |                                                                                                                                                                                                                                                                                                                 |
| Desorption | Instruments                | Turbomatrix, PerkinElmer, USA                                                                                 | Thermal desorption unit (TDU) and cooled injection system (CIS) (Gerstel, Germany)                                                                                                                                                                                                                              |
|            | Steps-Temperature-duration | 250 °C for 15 min                                                                                             | TDU: 220 °C for 5 min<br>CIS: −10 °C for 1 min ramped to 270 °C at a heating rate of 12 °C per second and held for 2.5 min in spitless mode                                                                                                                                                                     |
| Analysis   | Instruments                | Clarus 680, Perkin Elmer, USA                                                                                 | Agilent 7890B gas chromatograph                                                                                                                                                                                                                                                                                 |
|            | Columns                    | Elite-5-ms (60 m × 0.25 mm × 1 µm, Perkin Elmer, USA)                                                         | ZB-5MSplus fused silica capillary column (30 m × 0.25 mm × 0.25 µm, Phenomenex, Torrance, CA, USA)                                                                                                                                                                                                              |
|            | Gas-flow rate              | Helium-constant pressure of 30 psi                                                                            | Helium-0.8 mL/min                                                                                                                                                                                                                                                                                               |
|            | Oven program               | 9 min at 40 °C, a ramp at 15 °C/min until 90 °C hold 4 min, then a ramp at 10 °C/min until 250 °C hold 15 min | 4 min at 40 °C, a ramp at 6 °C/min until 300 °C, hold 1 min                                                                                                                                                                                                                                                     |
| Detection  | Instruments                | SQ8T model mass spectrometer (Perkin Elmer, USA)                                                              | 5977A mass spectrometer (Agilent, Santa Clara, CA, USA)                                                                                                                                                                                                                                                         |
|            | Mode                       | Electronic impact at 70 eV                                                                                    | Electronic impact at 70 eV                                                                                                                                                                                                                                                                                      |
|            | Acquired mass              | 20 to 350 amu                                                                                                 | 33 to 350 amu                                                                                                                                                                                                                                                                                                   |
|            | Analysis type              | Scan mode                                                                                                     | Scan mode                                                                                                                                                                                                                                                                                                       |
|            | Identification             | NIST library                                                                                                  | NIST library                                                                                                                                                                                                                                                                                                    |

<sup>a</sup> SBSE = Stir Bar Sorptive Extraction, amu = atomic mass unit. <sup>\*</sup> Stir bars were conditioned prior to use according to the manufacturer's instructions.

The generation of the controlled atmosphere is conducted as follows: volatilization of liquid solution of 6 compounds (Table 3) in a 20 L/min air flow; then dilution of this air flow in a 417 L/min air flow to reach the required concentration (417 L/min is the maximum air flow that can be set in the dilution pilot). We use a high flow rate because of the Coriolis sampling constraints. Three samples were realized by each technique in duplicate experiment.

**Table 3.** Information about substances contained in the solution used for air-controlled sampling.

| CAS       | Substances        | Family   | Molecular Weight (g/mol) | Concentration (mg/m <sup>3</sup> ) | References <sup>1</sup> |
|-----------|-------------------|----------|--------------------------|------------------------------------|-------------------------|
| 78-93-3   | 2-butanone        | Ketone   | 72.11                    | 3.59                               | VWR (25643.294)         |
| 108-88-3  | Toluene           | Aromatic | 92.14                    | 1.95                               | Supelco (1.00849.2500)  |
| 124-13-0  | Octanal           | Aldehyde | 128.21                   | 2.76                               | Acros (199481000)       |
| 104-76-7  | 2-ethyl-1-hexanol | Alcohol  | 130.23                   | 2.81                               | Acros (118530010)       |
| 111-87-5  | 1-octanol         | Alcohol  | 130.23                   | 3.69                               | Fluka (74852)           |
| 1120-21-4 | Undecane          | Alkane   | 156.31                   | 2.48                               | Acros (140665000)       |

<sup>1</sup> All the standards have a purity higher than 99% (chromatography grade).

Sampling using the Coriolis micro and Nalophan<sup>®</sup> bags were performed as described previously. Each sampling was performed in triplicate. The volumes of liquid that were withdrawn with the Coriolis micro were 15, 22, and 30 mL, respectively.

#### 4.3. Estimated Data

The physico-chemical properties of each compound were determined using the Estimation Programs Interface Suite 4.11 (EPI Suite<sup>TM</sup>) [56].

KOWINTM was used for estimating log  $K_{ow}$  with an atom/fragment contribution method, HENRYWINTM was utilized to determine Henry's constant using the group and bound contribution. The vapor pressure and boiling points for each compound has been determined by the MPBPVPTM tool and KOAWINTM was used to estimate the  $K_{oa}$ , using the ration of the log  $K_{ow}$  from KOWINTM and the dimensionless Henry's Law constant from HENRYWINTM.

To obtain all these values, we entered the Simplified Molecular-Input Line-Entry System (SMILES) of each molecule in the software.

#### 4.4. Statistical Analysis

A statistical analysis and drawing were performed using the R software, version 4.1.1 by R core team [57]. The data were analyzed by the Shapiro-Wilk test to determine if our data followed a normal distribution. If data followed a normal distribution a *t*-test of Student for unpaired data was applied, if not a Test U from Mann-Whitney was chosen. For all statistical tests a P value threshold of 0.05 was used. A Venn diagram was realized using the software interactivenn.net (date of access: 01/06/2022) [58].

**Supplementary Materials:** The following supporting information can be downloaded at: <https://www.mdpi.com/article/10.3390/molecules27196462/s1>, Table S1: physico-chemicals properties of molecules found with the Coriolis micro; Table S2: physico-chemicals properties of molecules found with Nalophan<sup>®</sup> bags. Table S3: details of statistical tests.

**Author Contributions:** Conceptualization, data curation, methodology, writing—original draft A.C. (Audrey Courtier); funding acquisition, project administration, supervision, A.C. (Axelle Cadere) and S.B.; validation, S.C.; writing—review and editing, B.R., S.C., A.C. (Axelle Cadere) and S.B. All authors have read and agreed to the published version of the manuscript.

**Funding:** This research was funded by the Occitanie region and European Regional Development Fund (FEDER).

**Acknowledgments:** The authors thank A. Comte for proofreading.

**Conflicts of Interest:** The authors declare no conflict of interest.

**Sample Availability:** Samples of the compounds are available from the authors.

## References

1. World Health Organization. Air Pollution. Available online: <https://www.who.int/health-topics/air-pollution> (accessed on 1 August 2022).
2. Khreis, H.; Cirach, M.; Mueller, N.; de Hoogh, K.; Hoek, G.; Nieuwenhuijsen, M.J.; Rojas-Rueda, D. Outdoor Air Pollution and the Burden of Childhood Asthma across Europe. *Eur. Respir. J.* **2019**, *54*, 1802194. [CrossRef] [PubMed]
3. Xu, X.; Nie, S.; Ding, H.; Hou, F.F. Environmental Pollution and Kidney Diseases. *Nat. Rev. Nephrol.* **2018**, *14*, 313–324. [CrossRef] [PubMed]
4. Turner, M.C.; Andersen, Z.J.; Baccarelli, A.; Diver, W.R.; Gapstur, S.M.; Pope, C.A.; Prada, D.; Samet, J.; Thurston, G.; Cohen, A. Outdoor Air Pollution and Cancer: An Overview of the Current Evidence and Public Health Recommendations. *CA Cancer J. Clin.* **2020**, *70*, 460–479. [CrossRef] [PubMed]
5. Kurt, O.K.; Zhang, J.; Pinkerton, K.E. Pulmonary Health Effects of Air Pollution. *Curr. Opin. Pulm. Med.* **2016**, *22*, 138–143. [CrossRef]
6. US EPA Organization. What Is the Definition of VOC? Available online: <https://www.epa.gov/air-emissions-inventories/what-definition-voc> (accessed on 9 November 2021).
7. Raffy, G.; Mercier, F.; Blanchard, O.; Derbez, M.; Dassonville, C.; Bonvallot, N.; Glorennec, P.; Le Bot, B. Semi-Volatile Organic Compounds in the Air and Dust of 30 French Schools: A Pilot Study. *Indoor Air* **2017**, *27*, 114–127. [CrossRef]
8. Jia, L.; Chu, J.; Li, J.; Ren, J.; Huang, P.; Li, D. Formaldehyde and VOC Emissions from Plywood Panels Bonded with Bio-Oil Phenolic Resins. *Environ. Pollut.* **2020**, *264*, 114819. [CrossRef]
9. Zhang, T.; Li, G.; Yu, Y.; Ji, Y.; An, T. Atmospheric Diffusion Profiles and Health Risks of Typical VOC: Numerical Modelling Study. *J. Clean. Prod.* **2020**, *275*, 122982. [CrossRef]
10. Sazakli, E.; Leotsinidis, M. Odor Nuisance and Health Risk Assessment of VOC Emissions from a Rendering Plant. *Air Qual. Atmos. Health* **2021**, *14*, 301–312. [CrossRef]
11. Mo, Z.; Lu, S.; Shao, M. Volatile Organic Compound (VOC) Emissions and Health Risk Assessment in Paint and Coatings Industry in the Yangtze River Delta, China. *Environ. Pollut.* **2021**, *269*, 115740. [CrossRef]
12. Jantunen, M.; Jaakkola, J.J.K.; Krzyzanowski, M. (Eds.) *Assessment of Exposure to Indoor Air Pollutants*; WHO Regional Publications; World Health Organization Regional Office for Europe: Copenhagen, Denmark, 1997.
13. US EPA. Technical Overview of Volatile Organic Compounds. Available online: <https://www.epa.gov/indoor-air-quality-iaq/technical-overview-volatile-organic-compounds#8> (accessed on 27 August 2021).
14. OMS. *Health Effects of Particulate Matter*; World Health Organization Regional Office for Europe: Copenhagen, Denmark, 2013.
15. Arh, G.; Klasinc, L.; Veber, M.; Pompe, M. Calibration of Mass Selective Detector in Non-Target Analysis of Volatile Organic Compounds in the Air. *J. Chromatogr. A* **2011**, *1218*, 1538–1543. [CrossRef]
16. Woolfenden, E. Sorbent-Based Sampling Methods for Volatile and Semi-Volatile Organic Compounds in Air. *J. Chromatogr. A* **2010**, *1217*, 2674–2684. [CrossRef]
17. Krupa, S.V.; Legge, A.H. Passive Sampling of Ambient, Gaseous Air Pollutants: An Assessment from an Ecological Perspective. *Environ. Pollut.* **2000**, *107*, 31–45. [CrossRef]
18. Seethapathy, S.; Górecki, T.; Li, X. Passive Sampling in Environmental Analysis. *J. Chromatogr. A* **2008**, *1184*, 234–253. [CrossRef]
19. Crump, D. Application of Diffusive Samplers. In *Organic Indoor Air Pollutants*; John Wiley & Sons, Ltd.: New York, NY, USA, 2009; pp. 47–63. [CrossRef]
20. Schieweck, A. Very Volatile Organic Compounds (VVO) as Emissions from Wooden Materials and in Indoor Air of New Prefabricated Wooden Houses. *Build. Environ.* **2021**, *190*, 107537. [CrossRef]
21. Woolfenden, E. Sorbent-Based Sampling Methods for Volatile and Semi-Volatile Organic Compounds in Air. Part 2. Sorbent Selection and Other Aspects of Optimizing Air Monitoring Methods. *J. Chromatogr. A* **2010**, *1217*, 2685–2694. [CrossRef]
22. McDonald, J.D.; Zielinska, B.; Fujita, E.M.; Sagebiel, J.C.; Chow, J.C.; Watson, J.G. Fine Particle and Gaseous Emission Rates from Residential Wood Combustion. *Environ. Sci. Technol.* **2000**, *34*, 2080–2091. [CrossRef]
23. Ibáñez, M.; Sancho, J.V.; Hernández, F.; McMillan, D.; Rao, R. Rapid Non-Target Screening of Organic Pollutants in Water by Ultrapformance Liquid Chromatography Coupled to Time-of-Light Mass Spectrometry. *TrAC Trends Anal. Chem.* **2008**, *27*, 481–489. [CrossRef]
24. Gosetti, F.; Mazzucco, E.; Gennaro, M.C.; Marengo, E. Contaminants in Water: Non-Target UHPLC/MS Analysis. *Environ. Chem. Lett.* **2016**, *14*, 51–65. [CrossRef]

25. Ruff, M.; Mueller, M.S.; Loos, M.; Singer, H.P. Quantitative Target and Systematic Non-Target Analysis of Polar Organic Micro-Pollutants along the River Rhine Using High-Resolution Mass-Spectrometry – Identification of Unknown Sources and Compounds. *Water Res.* **2015**, *87*, 145–154. [CrossRef]
26. Picó, Y.; El-Sheikh, M.A.; Alfarhan, A.H.; Barceló, D. Target vs Non-Target Analysis to Determine Pesticide Residues in Fruits from Saudi Arabia and Influence in Potential Risk Associated with Exposure. *Food Chem. Toxicol.* **2018**, *111*, 53–63. [CrossRef]
27. Angeles, L.F.; Halwatura, L.M.; Antle, J.P.; Simpson, S.; Jaraula, C.M.B.; Aga, D.S. In Silico Supported Nontarget Analysis of Contaminants of Emerging Concern: Increasing Confidence in Unknown Identification in Wastewater and Surface Waters. *ACS EST Water* **2021**, *1*, 1765–1775. [CrossRef]
28. Pourchet, M.; Debrauwer, L.; Klanova, J.; Price, E.J.; Covaci, A.; Caballero-Casero, N.; Oberacher, H.; Lamoree, M.; Damont, A.; Fenaille, F.; et al. Suspect and Non-Targeted Screening of Chemicals of Emerging Concern for Human Biomonitoring, Environmental Health Studies and Support to Risk Assessment: From Promises to Challenges and Harmonisation Issues. *Environ. Int.* **2020**, *139*, 105545. [CrossRef]
29. Wilczyńska-piliszek, A.J.; Piliszek, S.; Falandysz, J. Estimation of  $K_{OA}$  Values of 209 Polychlorinated *Trans*-Azobenzenes by PM6 and DFT Methods. *J. Environ. Sci. Health Part B* **2012**, *47*, 562–570. [CrossRef]
30. Xu, S.; Kropscott, B. Octanol/Air Partition Coefficients of Volatile Methylsiloxanes and Their Temperature Dependence. *J. Chem. Eng. Data* **2013**, *58*, 136–142. [CrossRef]
31. Collège National d'Experts en Environnement de l'Industrie Chimique. D'Experts en E. de l'Industrie C. Les Propriétés Environnementales Des Substances. 2010. Available online: <https://docplayer.fr/56425504-Les-proprietes-environnementales-des-substances.html> (accessed on 6 January 2022).
32. Tissier, C.; Morvan, C.; Bocquené, G. *Les Substances Prioritaires De La Directive Cadre Sur L'eau (DCE)*; Institut Ifremer: Brest, France, 2005; Volume 92.
33. Shaw, P.E. Review of Quantitative Analyses of Citrus Essential Oils. *J. Agric. Food Chem.* **1979**, *27*, 246–257. [CrossRef]
34. Hausch, B.J.; Lorjaroenphon, Y.; Cadwallader, K.R. Flavor Chemistry of Lemon-Lime Carbonated Beverages. *J. Agric. Food Chem.* **2015**, *63*, 112–119. [CrossRef]
35. Jang, M.; Kamens, R.M. Atmospheric Secondary Aerosol Formation by Heterogeneous Reactions of Aldehydes in the Presence of a Sulfuric Acid Aerosol Catalyst. *Environ. Sci. Technol.* **2001**, *35*, 4758–4766. [CrossRef]
36. Lee, A.K.Y.; Li, Y.J.; Lau, A.P.S.; Chan, C.K. A Re-Evaluation on the Atmospheric Significance of Octanal Vapor Uptake by Acidic Particles: Roles of Particle Acidity and Gas-Phase Octanal Concentration. *Aerosol Sci. Technol.* **2008**, *42*, 992–1000. [CrossRef]
37. International Agency for Research on Cancer (IARC). Pyridine (IARC Summary & Evaluation, Volume 77, 2000). Available online: <https://inchem.org/documents/iarc/vol77/77-16.html> (accessed on 6 January 2022).
38. Ratiu, I.A.; Al-Suod, H.; Ligor, M.; Monedeiro, F.; Buszewski, B. Effects of Growth Conditions and Cultivability on the Content of Cyclitols in *Medicago Sativa*. *Int. J. Environ. Sci. Technol.* **2021**, *18*, 33–48. [CrossRef]
39. Zalut, O.A.; Elsayed, M.A. A Study on Microwave Removal of Pyridine from Wastewater. *J. Environ. Chem. Eng.* **2013**, *1*, 137–143. [CrossRef]
40. Lataye, D.H.; Mishra, I.M.; Mall, I.D. Removal of Pyridine from Aqueous Solution by Adsorption on Bagasse Fly Ash. *Ind. Eng. Chem. Res.* **2006**, *45*, 3934–3943. [CrossRef]
41. USEPA. List of Lists: Consolidated List of Chemicals Subject to the Emergency Planning and Community Right-To-Know Act (EPCRA), Comprehensive Environmental Response, Compensation and Liability Act (CERCLA) and Section 112(r) of the Clean Air Act. 2021, 131. Available online: <https://www.epa.gov/epcra/consolidated-list-lists-under-epcracerclaa-ss112r-april-2022-version> (accessed on 6 January 2022).
42. Wang, X.; Wei, Y.; Zhang, H.; Bao, L.; He, M.; Yuan, S. Understanding the Properties of Methyl Vinyl Ketone and Methacrolein at the Air-Water Interface: Adsorption, Heterogeneous Reaction and Environmental Impact Analysis. *Chemosphere* **2021**, *283*, 131183. [CrossRef] [PubMed]
43. Biesenthal, T.A.; Shepson, P.B. Observations of Anthropogenic Inputs of the Isoprene Oxidation Products Methyl Vinyl Ketone and Methacrolein to the Atmosphere. *Geophys. Res. Lett.* **1997**, *24*, 1375–1378. [CrossRef]
44. ECHA. Substance Information—ECHA. Available online: <https://echa.europa.eu/fr/substance-information/-/substanceinfo/100.001.055> (accessed on 6 January 2022).
45. Ramírez, N.; Cuadras, A.; Rovira, E.; Borrull, F.; Marcé, R.M. Chronic Risk Assessment of Exposure to Volatile Organic Compounds in the Atmosphere near the Largest Mediterranean Industrial Site. *Environ. Int.* **2012**, *39*, 200–209. [CrossRef]
46. Li, K.; Liang, T.; Wang, L. Risk Assessment of Atmospheric Heavy Metals Exposure in Baotou, a Typical Industrial City in Northern China. *Environ. Geochem. Health* **2016**, *38*, 843–853. [CrossRef]
47. Hu, Y.-J.; Bao, L.-J.; Huang, C.-L.; Li, S.-M.; Zeng, E.Y. A Comprehensive Risk Assessment of Human Inhalation Exposure to Atmospheric Halogenated Flame Retardants and Organophosphate Esters in an Urban Zone. *Environ. Pollut.* **2019**, *252*, 1902–1909. [CrossRef]
48. Gholizadeh, A.; Taghavi, M.; Moslem, A.; Neshat, A.A.; Lari Najafi, M.; Alahabadi, A.; Ahmadi, E.; Ebrahimi aval, H.; Asour, A.A.; Rezaei, H.; et al. Ecological and Health Risk Assessment of Exposure to Atmospheric Heavy Metals. *Ecotoxicol. Environ. Saf.* **2019**, *184*, 109622. [CrossRef]

49. Guida, Y.; Pozo, K.; de Carvalho, G.O.; Capella, R.; Targino, A.C.; Torres, J.P.M.; Meire, R.O. Occurrence of Pyrethroids in the Atmosphere of Urban Areas of Southeastern Brazil: Inhalation Exposure and Health Risk Assessment. *Environ. Pollut.* **2021**, *290*, 118020. [CrossRef]
50. Skjevrak, I.; Brede, C.; Steffensen, I.-L.; Mikalsen, A.; Alexander, J.; Fjeldal, P.; Herikstad, H. Non-Targeted Multi-Component Analytical Surveillance of Plastic Food Contact Materials: Identification of Substances Not Included in EU Positive Lists and Their Risk Assessment. *Food Addit. Contam.* **2005**, *22*, 1012–1022. [CrossRef]
51. Hu, Y.; Du, Z.; Sun, X.; Ma, X.; Song, J.; Sui, H.; Debrah, A.A. Non-Targeted Analysis and Risk Assessment of Non-Volatile Compounds in Polyamide Food Contact Materials. *Food Chem.* **2021**, *345*, 128625. [CrossRef]
52. McCord, J.P.; Groff, L.C.; Sobus, J.R. Quantitative Non-Targeted Analysis: Bridging the Gap between Contaminant Discovery and Risk Characterization. *Environ. Int.* **2022**, *158*, 107011. [CrossRef]
53. SMBT. *SAGE Des Bassins Versants de La Lagune de Thau et de l'étang d'Ingril*; SMBT: Sète, France, 2016.
54. DREAL. Réseau de Suivi de la Lagune de Thau—DREAL Occitanie. Available online: <https://www.occitanie.developpement-durable.gouv.fr/reseau-de-suivi-de-la-lagune-de-thau-r1326.html> (accessed on 13 June 2022).
55. Berrou, K.; Dunyach-Remy, C.; Lavigne, J.-P.; Roig, B.; Cadiere, A. Comparison of Stir Bar Sorptive Extraction and Solid Phase Microextraction of Volatile and Semi-Volatile Metabolite Profile of *Staphylococcus Aureus*. *Molecules* **2020**, *25*, 55. [CrossRef]
56. US EPA. *Estimation Programs Interface Suite™ for Microsoft®Windows*; version 4.11; United States Environmental Protection Agency: Washington, DC, USA, 2021.
57. Bunn, A.; Korpela, M. Crossdating in DplR. 2014, 12. Available online: <https://mran.microsoft.com/snapshot/2019-04-29/web/packages/dplR/vignettes/xdate-dplR.pdf> (accessed on 6 January 2022).
58. Heberle, H.; Meirelles, G.V.; da Silva, F.R.; Telles, G.P.; Minghim, R. InteractiVenn: A Web-Based Tool for the Analysis of Sets through Venn Diagrams. *BMC Bioinform.* **2015**, *16*, 169. [CrossRef]